



Grace Care series

D2w

TWO FUNCTIONS MANUAL BED

Registration code: SK-CD1




SAIKANG/Medical

Healthier
your life

D2w

FEATURES AND HIGHLIGHTS

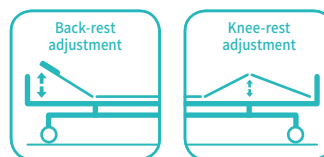

Technical parameters

▪ External size (LxWxH)	2245x1000x510 mm
▪ Mattress platform	1925x830 mm
▪ Safe working load	250 kg
▪ Back-rest adjustment	0-75° (±10°)
▪ Knee-rest adjustment	0-35° (±10°)



Technical configuration

☒ Manual Cranks	2pcs
☒ 5" Hardshell Castors	4pcs
☒ Foldable Side Rail with ProErgo-Handle™	1set
☒ ABS Grace Bed Ends with Safe Lock	1set
☒ IV Pole	1set
☒ IV Pole Prevision	4pcs
☒ Drainage Hook	2pcs
☒ Shoes Holder	1set



Optional

☐ Monkey Pole	
☐ Balkan Traction Frame	
☐ Directional Wheel	



Design

Extra long and height head foot board, with curve surrounded design, provide comforting feeling for patients.

**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa Virgin ABS material.

**Steels | Painting**

Baosteel from Fortune Global 500 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth wedding Solid bed base made of 5.0 mm stamped rotating plate.

**Hand Control**

Retractable hand crank, 100,000 times remaining intact
Extreme protection device, sounds when tipped to the max position,
Strong pallets to increase the bed height adjustment.

**Platform**

Stamping forming platform, smooth surface no welding scar, bearing 250KG load. Option resin X-ray translucent platform.

**Guardrail**

L-shaped foldable guardrail (Patent) with ProErgo-Handle™, patented parts. with red and green mark to indicated if the siderail is locked or not. Bended tube, solid installing, 0 shaking, anti-pinch, protect patient's safety.

**Head & Foot Board**

Head and foot board can be reversible, extra height design, size: 1000*563mm; ABS Grace Bed Ends with Safe Lock, automatic lock without any manual operation on lock, on/off removable mark on hook lock.

**Castors**

5" hardshell casters
TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts.
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 time.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic hook.

**Shoes Holder**

Shoes holder: 800*400mm, adjustable





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎ +86 18021231901 · 🕒 24/7
✉ sales@chinasaikang.com
🌐 www.saikangmedical.com



Company video



General catalogue