



GT3k

MANUAL ORTHOPEDIC TRACTION BED

Registration code: SKD-4C



SAIKANG Medical

Healthier
your life

GT3k

FEATURES AND HIGHLIGHTS



Technical parameters

▪ External size (LxWxH)	2270x970x500 mm
▪ Mattress platform	1925x830 mm
▪ Safe working load	250 kg
▪ Back-rest adjustment	0-75° ($\pm 10^\circ$)
▪ Knee-rest adjustment	0-35° ($\pm 10^\circ$)
▪ Trendelenburg	0-6°



Technical configuration

☒ Manual cranks	3 pcs
☒ 5" double side castors, central locking	4 pcs
☒ Central locking pedal	2 pcs
☒ Aluminium alloy ABS head and foot board with safe lock	1 set
☒ Aluminium alloy foldable side rail	1 set
☒ IV Pole	1 set
☒ IV Pole Prevision	4 pcs
☒ Drainage Hook	2 pcs
☒ Aluminium balkan frame	1 set



Steels | Painting

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding Solid bed base made of 5.0 mm stamped rotating plate.

**Balkan Frame**

Aluminum Balkan Frame, which can adjust the longitudinal, lateral and height positions in all directions, and can accurately adjust the most comfortable posture of the patient.

**Hand Control**

Retractable hand crank, 100,000 times remaining intact
Extreme protection device, sounds when tipped to the max position, Strong pallets to increase the bed.

**Bumper Cover**

Reduce the impact of the collision, safe and slinky, protect the bed and patient during transportation.

**Platform**

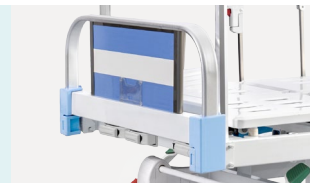
Stamping forming platform, smooth surface no welding scar, bearing 250KG load. Option resin X-ray translucent platform.

**Guardrail**

Aluminum alloy side rail, 6-staff, collapsed type with safe lock.

**Head & foot board**

Aluminium alloy flat tube, MDF material fire-proof plate. size of the head&foot board can be customized.

**Wheels**

Castors with TPR tire, no worn-out after running 30KM, united forming without bolts.
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.

**Brake**

Three position braking system, transporting can maintain a straight line.
ABS cover on castors joint, prevent the dust or hair stuck in. Easy to clean.

**Anti-slip bed feet cover**

Non-skid, more stable in the ground.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic hook.





CONSTANT
INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE
BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎ +86 18021231901 · 🕒 24/7
✉ sales@chinasaikang.com
🌐 www.saikangmedical.com



Company video



General catalogue