

SKB041-1


PATIENT TRANSPORTATION TROLLEY

Registration code: SKM-F



SKB041-1

FEATURES AND HIGHLIGHTS

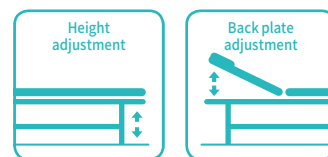

Technical parameters

▪ External size (LxW):	1900x640 mm
▪ Table height adjustment:	550-850 mm
▪ Back plate adjustment:	0-70° (±10°)

Technical configuration

<input checked="" type="checkbox"/> Crank system	1 set
<input checked="" type="checkbox"/> Back plate air spring	1 set
<input checked="" type="checkbox"/> ABS side guardrail	1 set
<input checked="" type="checkbox"/> 6" castors with central locking system	4 pcs
<input checked="" type="checkbox"/> Directional wheel	1 pc
<input checked="" type="checkbox"/> Central locking pedal	2 pcs
<input checked="" type="checkbox"/> Directional pedal	2 pcs
<input checked="" type="checkbox"/> IV pole	1 set
<input checked="" type="checkbox"/> IV pole prevision	2 pcs
<input checked="" type="checkbox"/> Oxygen bottle holder	1 pc
<input checked="" type="checkbox"/> Transfer mattress	1 pc
<input checked="" type="checkbox"/> Strips	2 pcs



Metal material

Baosteel from Fortune Global 500
Panasonic Robotic ensure 360° full
smooth welding.

**Bed frame&chassis**

Both with 30*50mm square tube,
more strong and stable.

**Painting**

11 process epoxy painting, ASTM
testing anti-bacterial, paint thickness
0.12 mm, brightness 60°. Painting can
resist 50 kg impact.

**Bed platform&guardrail**

Bed platform and guardrail are made
of high quality PP material.

**Plastic**

Environmental plastic can be
sterilized in 100°C, tensile upto 30
MPa.

**Back lift**

Back lift can be up to 70 degrees,
make patient comfortable when
needed.

**Double side castors**

Four double side castors of 6 inch
diameter, simultaneously lockable by
the step pedal at both sides.
TPR tire no worn out after running
30 km, save anti-winding hard shell,
united forming without bolts.

**Directional center wheel**

Can be engaged to facilitate direction
control, enabling users to navigate
the stretcher easily in hectic and
dynamic environments.

**Telescopic IV pole**

Plastic combined with stainless
steel 4 hooks IV pole (Φ16-19), two
IV pole previsions seperated on the
head&foot part.

**Oxygen bottle holder**

Two positions for 2L, Φ105mm oxygen
bottle, one is under head board,
another is on the chassis with belt.

**Stainless steel handle**

Stainless steel handle, more durable.

**Strips**

With Strips to prevent patients from
possibly sliding off the stretcher.
Standard thin mattress, 2cm thick,
density: 28kg/m3, color: dark blue.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service

WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa


Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
 +86 18021231901 ·  24/7
 sales@chinasaikang.com
 www.saikangmedical.com



Company video



General catalogue